

MiniFinder® Veo

Product manual



Innehållsförteckning

MiniFinder Veo

| | |
|--|---|
| <u>Package Contents</u> | 1 |
| <u>Overview</u> | 1 |
| <u>Management</u> | 2 |
| <u>Getting Started</u> | 2 |
| <u>Explanation of Voice Messages</u> | 3 |
| <u>Functions & settings</u> | 4 |
| <u>Activate SOS Alarm</u> | 5 |
| <u>Cancel SOS Alarm</u> | 5 |
| <u>Activate No Motion Alarm</u> | 6 |
| <u>Cancel No Motion Alarm</u> | 6 |
| <u>Fall Alarm</u> | 7 |
| <u>Specifications</u> | 8 |

Package Contents

1 MiniFinder Veo

2 Charging pad

3 USB Charging Cable


4 Wall Charger



Overview

MiniFinder Veo is a modern personal safety alarm designed to increase safety and independence.

Please read this guide before you start using the device.

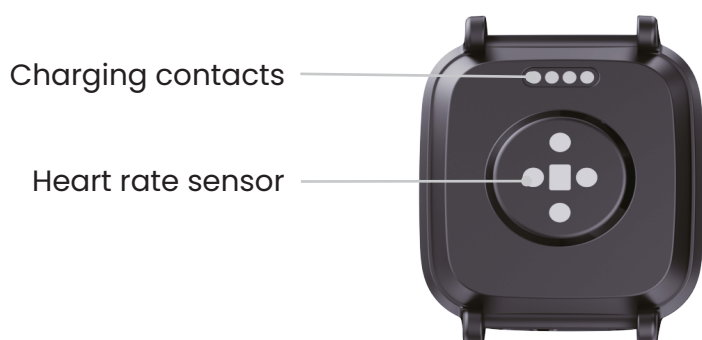

Management

1 Charge Before First Use

Charge the device for at least 3 hours before using it for the first time.

2 Charge the Device

To charge the device, place the magnetic charging pad on the back of the unit.

3 While Charging

During charging, the current battery percentage is displayed on the screen.

4 Power On the Device

To turn on the device, press and hold the front button until the MiniFinder logo appears.


5 Power Off the Device

To turn off the device, go to Settings > System > Power Off.

Note that the device cannot be turned off while it is charging.

Getting Started

To set up your MiniFinder Veo, start by charging the device. Connect the USB cable from the charging dock to the power adapter and plug it into a wall outlet.

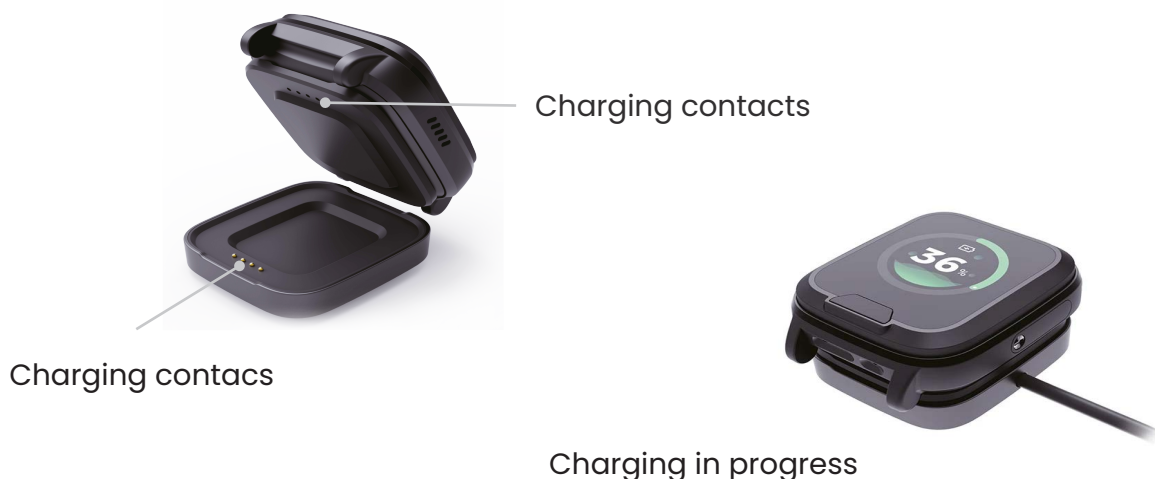

Explanation of Voice Messages

MiniFinder Veo can provide voice messages with information, such as charging status or different operating modes. The function is disabled by default, but the speaker volume and voice messages can be activated or adjusted in the app.

- A The battery in your MiniFinder is charging.**
Your MiniFinder Veo is connected to a power source and is currently charging.
- B The battery in your MiniFinder is low. Connect the charger**
The battery level of your MiniFinder Veo is low and should be charged. The message will repeat three times.
- C Your SOS alarm has been activated. Press the alarm button to cancel.**
The SOS alarm has been triggered. Press the alarm button to cancel the alarm before the call starts.
- D Your fall alarm has been activated. Press the alarm button to cancel.**
The fall alarm on your MiniFinder Veo has been activated. Press the alarm button to cancel the alarm.
- E Calling contact number...**
The device calls the preset emergency contacts in the order they are listed in the app. If the first emergency contact does not answer, the device will automatically proceed to the next contact in the list.

Functions & Settings

When MiniFinder Veo is turned on, swipe left or right to navigate through all available functions, such as the heart rate monitor, settings, and more.



Activate SOS Alarm

Press and hold the alarm button for 3 seconds. The device will vibrate and display the SOS alarm status on the screen. The device will then play a voice message::

“Your alarm has been activated”

The message is repeated three times. If the alarm is not canceled, the device will start calling the first preset contact.

Cancel SOS Alarm

The alarm can be canceled while the voice message is being played by pressing Cancel. If the alarm is not canceled, the device will continue to call the preset emergency contacts in the specified order. Once the calling sequence has started, the alarm cannot be canceled by default.



Activate No Motion Alarm

The No Motion alarm detects if the device has not registered any movement for an extended period. If no movement is detected during the set time, the No Motion alarm is activated. The device will then play a voice message:

“No movement detected, alarm will be activated”

The message is repeated three times. If the alarm is not canceled, the device will start calling the first preset contact.

The No Motion alarm is always active and will trigger if the device remains still for a preset period of time. The function is designed to detect extended periods of inactivity.

Activate No Motion Alarm

While the voice message is being played, the user can cancel the alarm by pressing the alarm button.

Once the device starts calling a contact, the alarm has been sent and can no longer be canceled.



Alarm Button

Fall Alarm

MiniFinder Veo can automatically detect a fall. When this occurs, the watch displays a warning message to the user.

While the voice message is being played, the user can cancel the alarm by pressing Cancel. If the alarm is not canceled within the preset time period, the device will call the emergency contacts.



Specifications

| | |
|--------------|-----------------------------|
| Dimensions | 47 mm x 42 mm x 15 mm |
| Weight | 38g |
| Battery Life | 48h with 3 minutes interval |
| Connectivity | LTE/4G (CatM1), GSM |
| SIM | Preinstalled |
| Voice Call | Yes |
| Fall sensor | Yes |
| Screen | 1.3" TFT (240 x 240) |
| GPS | GPS, GLONASS, Galileo |
| Waterproof | Yes, IPX7 |
| Button | SOS |

Memory

The device is equipped with a built-in 64 MB flash memory. This memory stores GPS data when GSM network coverage is unavailable (e.g., in rural areas, indoors, underground locations, etc.). The stored location data is automatically transmitted as soon as GSM coverage is restored.

Safety Declaration

This equipment has been tested and found to comply with safety certification requirements according to the specifications of EN Standard: EN60950-1: 2006 + A12: 2011.

Important Safety Instructions

Read these instructions carefully.

Keep these instructions for future reference.

Heed all warnings.

Follow all instructions.

Do not tamper with your MiniFinder®.

Do not use abrasive cleaners to clean your MiniFinder®.

Do not place your MiniFinder® in a dishwasher, washing machine, or dryer.

MiniFinder® should be used within -20°C to +60°C. Outside this range, the device is negatively affected

Do not use your MiniFinder® in a sauna or steam room.

Do not leave your MiniFinder® in direct sunlight for extended periods.

Do not place your MiniFinder® near an open flame.

Do not disassemble your MiniFinder® – it contains no user-serviceable parts.

The warranty will be void if the device is disassembled..

Recycling Information



The symbol on the product or its packaging indicates that this product must be disposed of separately from household waste at the end of its service life. Please be aware that it is your responsibility to dispose of electronic equipment properly at designated recycling centers, in order to help conserve natural resources.