

MiniFinder® Pico

Product Manual

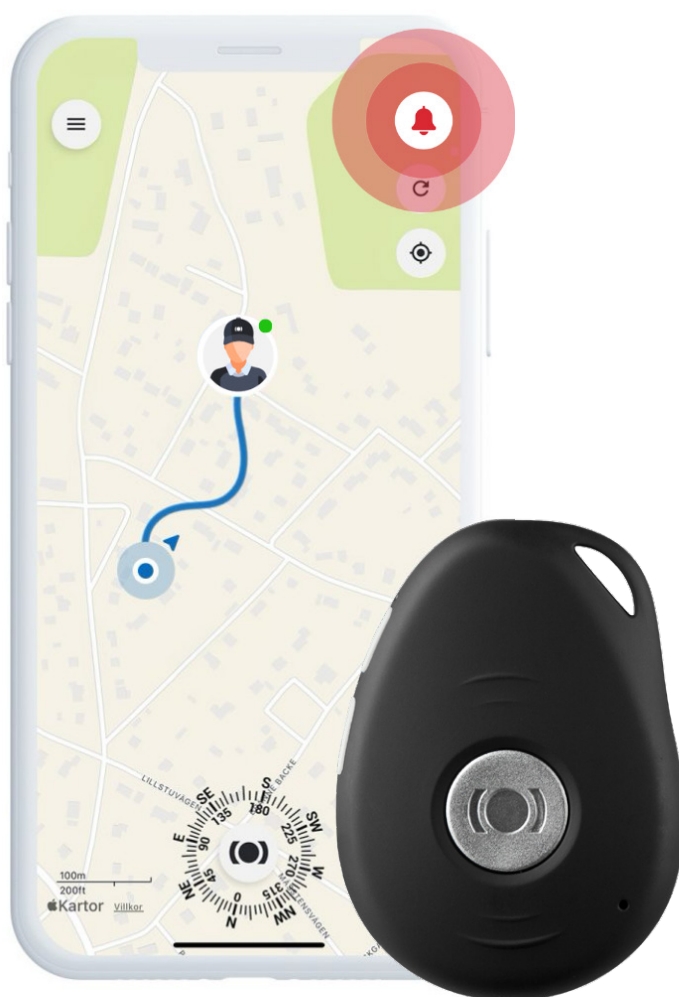


Table of contents

MiniFinder Pico

Package contents 1

Product Management 1

Explanation of LED Lights 2

Getting Started 2

Explanation of voice notifications 3

Activate SOS-Alarm 4

Cancel SOS-Alarm 4

Fall Alarm 5

Specifications 6

Package Contents




Product Management

- 1 Charge Before First Use**
Please charge the device for at least 4 hours before using it for the first time
- 2 While Charging**
The blue LED light will remain solid during charging.
Once the device is fully charged, the blue LED will turn off
- 3 Charging the Device**
Attach the magnetic charging cable to the charging pins on the device
- 4 Power On**
Press and hold the On/Off button for 1 second to turn on the device
- 5 Power Off**
Press and hold both SOS button and On/Off button for 2 seconds



Explanation of LED Lights

GREEN LED:

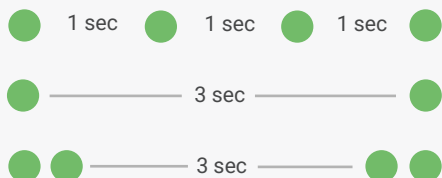

Rapid flash = Calling

Single flash = Searching for mobile network

Double flash = Connected to mobile network


RED LED:



Slow flash = Charging

Steady light = Fully charged

BLUE LED:

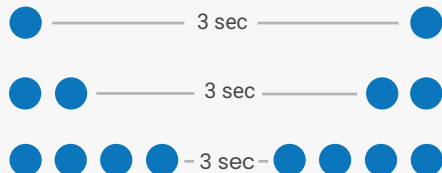

Slow flash = Searching for GPS signal

Double flash = GPS signal acquired

Blinks four times = Connected to WiFi

Getting Started – Download the MiniFinder

To begin using your **MiniFinder Pico**, you first need to download the **MiniFinder Live App**.

If you're viewing this guide on your mobile device, simply click the Google Play or App Store icons to go directly to the download page. Alternatively, open the App Store or Google Play, and search for "MiniFinder Live" to install the app manually.



Explanation of Voice Notifications

A Your mobile alarm is charging

Your MiniFinder Pico has been plugged in to a power source and is now charging

B Battery is low, please recharge your battery

Your MiniFinder Pico is running out of battery and you should charge it soon. This will be repeated 3 times

C Your alarm has been activated, press the SOS button to stop

The SOS-alarm has been activated. Press the MiniFinder button to cancel the alarm before it starts calling

D Fall alarm has been activated, press the SOS button to stop

Your MiniFinder Pico has activated the fall alarm. Press the MiniFinder button to cancel the alarm

E Calling contact number X

The device is now calling the authorized phone numbers you set up in the app. It will proceed with the next phone number if the current one doesn't respond.

Activate SOS Alarm

Press and hold the alarm button for 3 seconds. The device will vibrate and the device will then play a voice message:

“Your alarm has been activated”

The message is repeated three times. If the alarm is not canceled, the device will start calling the first preset contact.

Cancel SOS Alarm

The alarm can be canceled while the voice message is being played by pressing the SOS button. If the alarm is not canceled, the device will continue to call the preset emergency contacts in the specified order. Once the calling sequence has started, the alarm cannot be canceled by default.



3 Seconds
Alarm starts

>>>



Press to
cancel the alarm

Fall Alarm

MiniFinder Pico can automatically detect a fall. When this occurs, the watch displays a warning message to the user.

While the voice message is being played, the user can cancel the alarm by pressing the SOS button. If the alarm is not canceled within the preset time period, the device will call the emergency contacts.



Vibrates during
automatically fall alarm



Press to cancel
fall alarm

Specifications

Dimensions	60 mm x 41 mm x 16 mm
Weight	40g
Battery Life	36-40h with 5 minutes interval
Connectivity	LTE: B1, B3, B7, B8, B20
Power Source	5V DC 1.0A
SIM	Preinstalled
Voice call	Yes
Fall sensor	Yes
Sensor	Accelerometer, 3-axis

Safety Declaration

This equipment has been tested and found to comply with safety certification requirements according to the specifications of EN Standard: EN60950-1: 2006 + A12: 2011.

Important Safety Instructions

Read these instructions carefully.

Keep these instructions for future reference.

Heed all warnings.

Follow all instructions.

Do not tamper with your MiniFinder®.

Do not use abrasive cleaners to clean your MiniFinder®.

Do not place your MiniFinder® in a dishwasher, washing machine, or dryer.

Do not expose your MiniFinder® to extremely high or low temperatures.

Do not use your MiniFinder® in a sauna or steam room.

Do not leave your MiniFinder® in direct sunlight for extended periods.

Do not place your MiniFinder® near an open flame.

Do not attempt to disassemble your MiniFinder®; it contains no user-serviceable parts.

Recycling Information

The symbol on the product or its packaging indicates that this product must be disposed of separately from household waste at the end of its service life. Please be aware that it is your responsibility to dispose of electronic equipment properly at designated recycling centers, in order to help conserve natural resources.